

MYS-5320NE

3 MOTOR ELECTRIC HOSPITAL BED



The logo for Meyosis, featuring a stylized blue heart rate line followed by the word "Meyosis" in a bold, blue, sans-serif font. A registered trademark symbol (®) is located at the top right of the word.

www.meyosis.com

MYS-5320NE

GENERAL SPECIFICATIONS

- Backrest, legrest and knee movements, height adjustment, vascular positions are controlled by 3 motors
- Positions are controlled by hand control (on head side)
- The dual auto-regression system, reduces the pressure on the abdomen and the friction on the pelvis area and manual CPR
- Removable Polypropylene (PP) plastic head and foot ends with colour panels
- Split PP plastic side rails with soft drop real-ise mechanism and angle indicator
- 4 independent and detachable ABS lying sections
- Antibacterial electrostatic powder coated steel frame, open chasis

- ABS plastic covered chasis
- Ø125 mm, 360° swivel wheels with cross brake system
- Accessory adaptors on the corners, 1 height adjustable IV Pole
- Urine bag holders and belt hooks on both sides
- Protective plastic bumpers on corners
- Alternative color options for side rails
- 26 Density foam mattress with antibacterial, waterproof and zippered cover

OPTIONAL FEATURES

- Nurse control panel
- Central brake system
- Manual CPR
- Battery
- Side control panel



Movable Backrest

The sliding backrest protects the patient and prevents uncomfortable forces during position adjustments.

Position Safety

Side rails provide support for patients. The head-end grip bars make bed maneuvering easier.

Collision Prevention

Wall bumper cones protect the bed and room furniture from damage.



Hand Remote Control

Positions can be adjusted electronically via the buttons on the hand remote control.

Manual CPR (optional)

Easily accessible mechanical CPR handle on both sides.

Hygienic Mattress Bases

4-piece easy-to-clean ABS plastic mattress surface.

Height Adjustment Range

The 41 cm – 76 cm height range helps prevent patient falls and ensures a comfortable working environment for healthcare staff.

Cross Brake

360° rotatable castors with individual brake.

OPERATION ANGLES



Vascular
0° - 30°



Knee Movement



Back Movement
0° - 70°




Leg Movement
0° - 30°



Height

DIMENSIONS

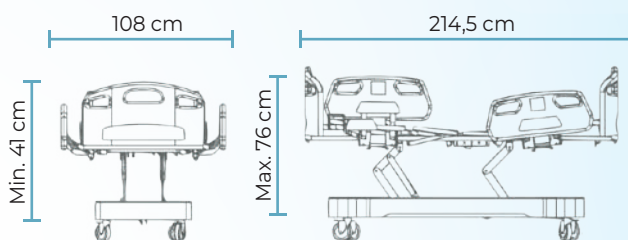

TECHNICAL SPECIFICATIONS



115 kg



Max. 250 kg



3 Motors



110-220V
50-60hz

